



Jackfish Electroless Deposition Kit

The information in this publication is provided for reference only. All information contained in this publication is believed to be correct and complete. PIKE Technologies, Inc. shall not be liable for errors contained herein nor for incidental or consequential damages in connection with the furnishing, performance, or use of this material. All product specifications, as well as the information contained in this publication, are subject to change without notice. This publication may contain or reference information and products protected by copyrights or patents and does not convey any license under the patent rights of PIKE Technologies, Inc. nor the rights of others. PIKE Technologies, Inc. does not assume any liability arising out of any infringements of patents or other rights of third parties.

This document contains confidential or proprietary information of PIKE Technologies, Inc. Neither this document nor the information herein is to be reproduced, distributed, used or disclosed, either in whole or in part, except as specifically authorized by PIKE Technologies, Inc.

PIKE Technologies, Inc. makes no warranty of any kind with regard to this material including, but not limited to, the implied warranties of merchantability and fitness for a particular purpose.

Copyright 1991-2025 by PIKE Technologies, Inc., Madison, WI 53719. Printed in the United States of America. All world rights reserved. No part of this publication may be stored in a retrieval system, transmitted, or reproduced in any way, including but not limited to, photocopy, photograph, magnetic or other record, without the prior written permission of PIKE Technologies, Inc.

Address comments to:

PIKE Technologies, Inc.
6125 Cottonwood Drive
Madison, WI 53719

Phone	(608) 274-2721
Fax	(608) 274-0103
Email	info@piketech.com
Web Site	www.piketech.com

March 7, 2025

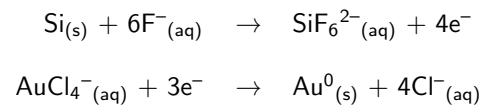
Contents

1	Introduction	1
2	The Electroless Deposition Kit	2
2.1	Packing List	2
3	Instructions	3
3.1	Deposition Solution Preparation	3
3.2	Crystal Preparation	3
3.3	Assembly of the Electroless Deposition Unit	3
3.4	Deposition	5
3.5	Fine-tuning the Process	5
4	Support	6
5	Also Available From Jackfish SEC	7

1 Introduction

The Jackfish SEC Electroless Deposition Fixture is designed to allow the end-user to chemically modify the reflecting plane of a face-angled crystal (FAC) with a metal film. Chemical-resistant Teflon and Viton materials prevent corrosive solution from causing damage to the sides of the crystal.

The electroless deposition process involves stripping the oxide layer on the Si surface followed by gold deposition via the dissolution of Si as follows:

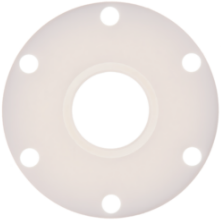



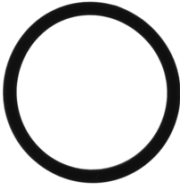



Assembled Electroless Deposition Kit.

2 The Electroless Deposition Kit

In order for you to quickly verify receipt of your accessory, we have included a packing list. Please inspect the package carefully. Contact PIKE Technologies immediately if any discrepancies are found.

2.1 Packing List

<p>Electroless deposition PTFE mask P.N. 162-4733 Qty. 1</p> 	<p>60 degree FAC holder, PEEK P.N. 162-4745 Qty. 1</p> 	<p>#8-32 screw, stainless steel, 1-1/2" long Qty. 8</p> 
<p>#8-32 nut, stainless steel Qty. 8</p> 	<p>-018 square profile FKM o-ring Qty. 2</p> 	<p>9/64" hex L-key Qty. 1</p> 

3 Instructions

3.1 Deposition Solution Preparation

CAUTION! The following solutions are extremely corrosive and present a significant health hazard. Use proper PPE and follow all relevant safety protocols when working with these chemicals.

The procedure described here was first described by Osawa in [1] and further modified in [2]. Recent use of this recipe on a FAC is described in [3].

Make all solutions with ultrapure water and analytical grade chemicals.

- 40% wt NH_4F – used to remove the oxide layer from the Si reflecting plane. Remake this solution regularly.
- Buffer Solution (0.3 M Na_2SO_4 , 0.1 $\text{Na}_2\text{S}_2\text{O}_3$, and 0.1 M NH_4Cl) – can be stored for several months of use.
- 2% HF – reduces the gold salt during the deposition process. Can be stored (in an appropriately compatible container) for several months of use.
- 0.015 - 0.030 M $\text{HAuCl}_4 \cdot 3\text{H}_2\text{O}$ – gold solution to be prepared fresh before deposition. Other gold salts may also be used.

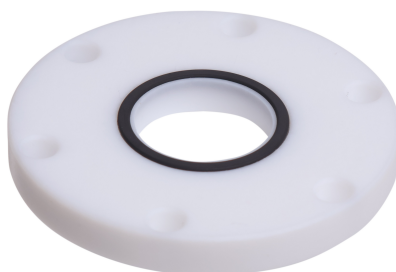
3.2 Crystal Preparation

Polishing the crystal is extremely important for high quality gold films. During the following steps, dust or debris on the crystal surface should be strictly avoided.


1. Using successively finer alumina slurries, polish the reflecting plane of the FAC for 20 minutes on each polish grade, rinsing with copious amounts of ultrapure water in between. *Polishing should be completed immediately prior to gold deposition!* We recommend two steps – one with 3 μm polish and another with 0.5 μm polish. Osawa and coworkers have reported using 0.03 μm polish (see [1]).
2. Sonicate the FAC in solvent (ethanol or acetone) and then in ultrapure water. Ensure that the FAC faces are carefully covered/protected during sonication to prevent damage to the crystal from contacting the walls of the sonication vessel.
3. Rinse with copious amounts of ultrapure water to remove all traces of solvent or debris. *Do not touch the reflecting plane during rinsing.*

3.3 Assembly of the Electroless Deposition Unit

Place the O-ring in the groove in the Teflon deposition mask.



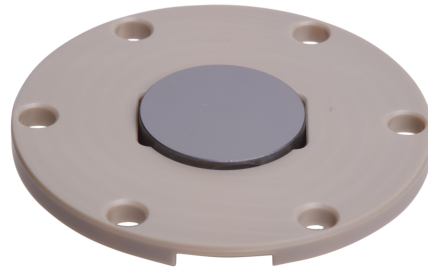Place the face-angled crystal (FAC) in the FAC holder and make sure the polished reflecting plane protrudes slightly above the top plane of the holder. Take care not to touch the reflecting plane.



Place the Teflon deposition mask, O-ring side down, onto the FAC. Ensure that the O-ring sits on top of the polished reflecting plane along the entire circumference of the FAC. Align the clearance holes on the FAC holder and deposition mask and place the screws in each of the six holes.



After placing a nut on each screw, hand tighten each of the six nuts. Then tighten with a screwdriver in a star pattern to ensure the deposition mask is even with the FAC holder. It is important to ensure that the mask does not shift from the FAC holder during this process.





Note: Once the assembly is complete, leak check the fixture by filling the well formed by the Teflon deposition mask and the surface of the FAC with water. Allow this to sit for at least 30 minutes and ensure there are no leaks prior to proceeding with corrosive solutions. This is especially important when using the Electroless Deposition Fixture for the first time.

The O-ring size was chosen to maximize the exposed area of the FAC reflecting plane. Any error in assembly can cause solution to leak onto the sides of the crystal which will have detrimental affects on surface quality.

When the assembly is complete, the long end of the screws should be coming out of the bottom of the FAC holder, as shown in the figure above. These will serve as a stand for the fixture.

Prepare a water bath held at 55°C with stirring. Ensure that the water level in the bath is sufficient to completely cover the bottom half of the FAC from the underside of the Electroless Deposition Fixture but low enough that it does not spill onto the FAC reflecting plane. Allow the FAC to come to the proper temperature before continuing.

3.4 Deposition

1. Freshly prepare the gold salt solution described in Section 3.1. We recommend about 1 mL of solution.
2. With a glass pipette, cover the FAC reflecting plane with NH_4F solution. Allow this solution to react for 3 minutes.
3. While waiting, add to a Teflon beaker (or other non-reactive beaker):
 - (a) Roughly 2 mL of the buffer solution
 - (b) 1.1 – 1.3 mL 2% HF solution. We recommend using a plastic syringe to measure this volume and mix it into the buffer solution.
4. When the three minutes have elapsed, add the gold salt solution to the beaker and mix well (this is the *plating solution*).
5. Empty the NH_4F from the Fixture.
6. Gently place 1.0 mL of plating solution on the FAC reflecting plane, avoiding bubbles.
7. Wait 8 to 10 seconds and gently remove 0.5 mL of plating solution from the FAC reflecting plane. Add another 0.5 mL.
8. Repeat Step 7 until all of the plating solution has been used.
9. Quench the reaction by removing the Electroless Deposition Fixture from the water bath and rinsing with at least 0.5 L of ultrapure water.
10. Disassemble the Electroless Deposition Fixture and rinse the gold-modified reflecting plane for at least 30 seconds with ultrapure water.
11. Carefully dispose of all waste and clean all components which have been exposed to corrosive solutions. We recommend prolonged immersion in ultrapure water in addition to any institutional protocols.

3.5 Fine-tuning the Process

This process requires significant practice and skill to achieve high-quality films. When troubleshooting, consider the cleanliness of the crystal and all chemicals and labware used in the process. If the film appears thick or uneven, then shorter deposition intervals in Step 7 of Section 3.4 can be tried, ranging from 4 to 6 seconds. Alternatively, the gold concentration can be lowered. If the film is still uneven, try a longer polishing procedure or with finer grade polish.

Gold films deposited in this fashion are extremely fragile – great care must be taken when assembling the Jackfish SEC J1 spectroelectrochemical cell to avoid damaging the gold layer.

4 Support

This manual and future updates are made available online at <https://jackfishsec.com/support>.

We have a blog which addresses some common questions such as preparing internal reflection elements for ATR-SEIRAS and selecting the optimal angle of incidence. Check it out at <https://jackfishsec.com/blog>.

We'd love to hear from you! Questions and feedback can be directed towards info@jackfishsec.com.

References

- (1) Miyake, H.; Ye, S.; Osawa, M. Electroless Deposition of Gold Thin Films on Silicon for Surface-Enhanced Infrared Spectroelectrochemistry. *Electrochemistry Communications* **2002**, *4* (12), 973–977, DOI: [10.1016/S1388-2481\(02\)00510-6](https://doi.org/10.1016/S1388-2481(02)00510-6).
- (2) Rosendahl, S. M.; Danger, B. R.; Vivek, J. P.; Burgess, I. J. Surface Enhanced Infrared Absorption Spectroscopy Studies of DMAP Adsorption on Gold Surfaces. *Langmuir* **2009**, *25* (4), 2241–2247, DOI: [10.1021/la803404u](https://doi.org/10.1021/la803404u).
- (3) Sigrist, J.; Lins, E.; Morhart, T.; Briggs, J.; Burgess, I. Optimization of a Commercial Variable Angle Accessory for Entry Level Users of Electrochemical Attenuated Total Reflection Surface Enhanced Infrared Absorption Spectroscopy (ATR-SEIRAS). *Applied Spectroscopy* **2019**, *73* (12), 1394–1402, DOI: [10.1177/0003702819858353](https://doi.org/10.1177/0003702819858353).

5 Also Available From Jackfish SEC

J1 Spectroelectrochemical Cell



Our flagship cell. Features a separate reference arm and glass construction for ultimate cleanliness.

J2 Spectroelectrochemical Cell



A completely sealed cell with a simple design. Compatible with both aqueous and volatile solvents.

J3 Easy Electrochemistry "EEL" Flow Cell



A cell with a small internal volume and flow capabilities, allowing researchers to investigate complex and dynamic surface phenomena.

J4 H-Cell



A cell with two glass compartments separated by a user-supplied membrane. Includes a quartz window for visual inspection of the working electrode.

Si Microgrooved Wafers



Uncoated microgrooved silicon wafers, available in a choice of two different face angles: 35° and 55° .



www.piketech.com | info@piketech.com | 608.274.2721 | Madison, WI USA

© Copyright PIKE Technologies, Inc. All world rights reserved.